



Taha Kimia Tajhiz Co.




Core Lab

Acoustic Velocity Core Holder- AVC Series Datasheet

Core Holder



Acoustic Velocity Core Holder- AVC Series

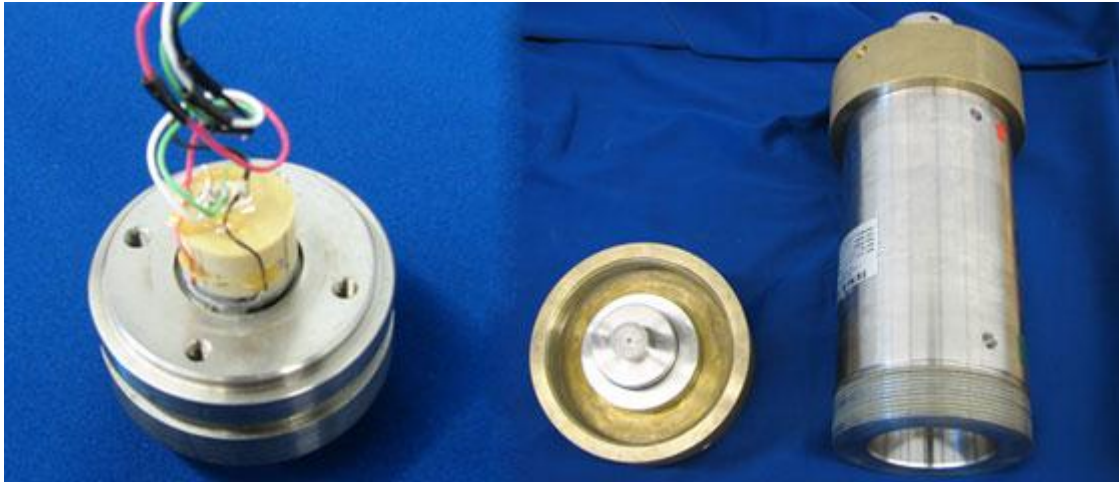

Acoustic velocity measurements are performed during logging measurements at reservoir locations. This core holder is used to test a core sample in the laboratory at reservoir conditions. The core holder contains a pair of insulated transducers.

The core sample is tri-axial loaded to insure proper surface contact between the sample and the transducers. The operator can flow fluids through the core sample in addition to performing the acoustic measurements. The operating pressure, temperature, and sample dimensions can be specified by the client.